

During the past three years our clients have received the benefits of our Esthetic Simplified Predictable removable system. ESP enables us to fabricate pre-fit prosthetics that minimize chairtime and ensure patient satisfaction. Join the "Removable Revolution", let ROE create your next denture case .

## ESP Complete Denture

## ESP Removables

### 1 Processed Base 1<sup>st</sup> Verification of Fit



For the bite appointment we deliver an Ivocap processed base enabling not just the pre-fit verification of the denture, but also providing a stable base for accurate bite registration. Any adjustments necessary to the final base should be made at this time.

### 2 Advanced Articulation & Set-up



Utilizing an average plane device, casts are mounted following anatomical landmarks. Our technicians utilize special tools to mechanically verify balanced lingual contact occlusion and curve of Spee and Wilson.

### 3 Try-In 2<sup>nd</sup> Verification of Fit




Now the patient can again verify the fit. Final processing involves only the acrylic around the teeth. Tissue bearing acrylic is untouched in the lab.

### 4 Beautiful, Patient Pre-Approved Finished Denture



Patient has pre-approved the final fit of the denture at both the bite and at try-in appointments. This will be a trouble-free appointment.

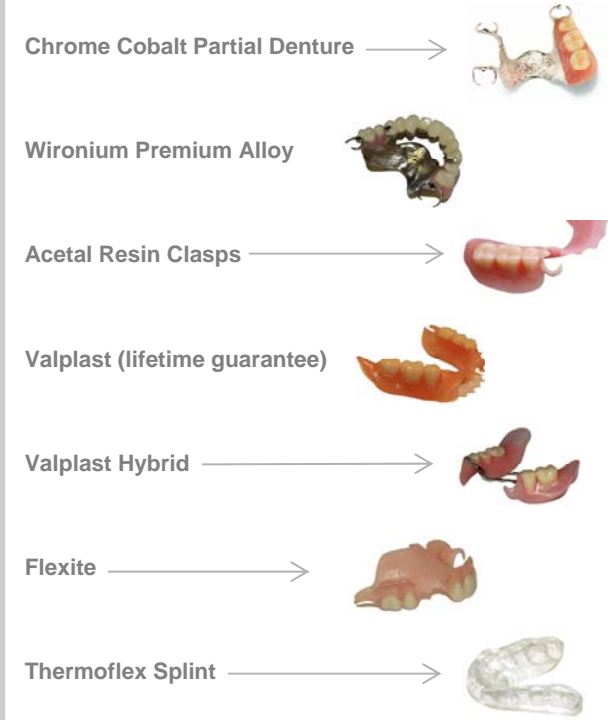

All removable appliances are fabricated on duplicate models. Partial dentures and splints are returned on the master cast for pre-seating verification of fit.